

# SALON LASER

## Scanning Hair Remover



This product **MUST** be used responsibly.  
Please read these Instructions carefully before use.



Refer to your Information Guide for:

- WARNINGS
- SECURITY BUTTON SEQUENCE
- HOW DO LASERS WORK?
- SUITABLE SKIN TYPES
- PATCH TEST
- HAIR GROWTH CYCLE
- WHY IS A SERIES OF TREATMENTS REQUIRED?
- YOUR QUESTIONS ANSWERED
- TROUBLESHOOTING

### INSTRUCTIONS

Please retain these instructions for future reference



808nm wavelength

Wearing laser safety goggles is unnecessary as the product has been designed with inbuilt safety features.

LASER SAFE  
TECHNOLOGY

PATENTS PENDING



CLASS 1 PRODUCT  
CONTAINING  
CLASS 3B LASER

I-N-LAS-UK/12


### warnings

Not suitable for Afro-Caribbean or black skin complexions.  
Not suitable for Grey Hair. See "Suitable Skin Types" in Information Guide.  
Laser hair removal should not cause discomfort. If you feel any discomfort or soreness the setting may be too high or your skin tone (colour) is not suitable for use with this product.  
Inbuilt safety features protect the eye being exposed to the laser, however, never aim the Wand at or near eyes.  
Do NOT use near flammable materials.  
Keep the product out of reach of children and pets.  
Reflection Hazard: Please remove all jewellery in the vicinity of the treatment area before use.  
When not in use, store the key responsibly in a safe place away from children.  
Do not treat the same hair repeatedly in any one session.  
Do not use after drinking alcohol or if you are feeling unwell or fatigued.  
Do not treat the following areas: tattoos, permanent make-up, pimples, warts, hair on black spots or moles, dark freckles, skin disease, birth marks, burns, lips, nipples, genital or anal areas, eczema, eyelashes, eyebrows, near the eyes, nasal or ear hair.  
This is a sensitive laser scanning product. Care should be taken to avoid damage to the mechanism in the Wand.

### cautions

Do NOT dismantle or repair. If you have any problems using the product, contact our Customer Care Department.  
Misuse of the unit can cause harm.  
Do NOT override any safety features or open the unit or wand.  
Only use this product with the Mains Adapter supplied.  
Keep Wand away from strong electromagnetic fields and magnets during use.


## components

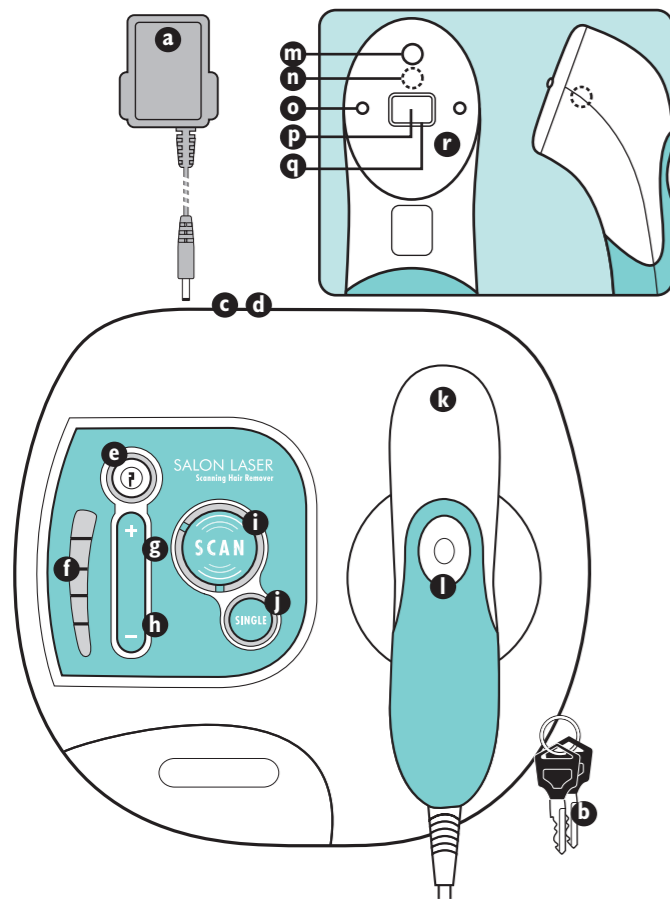

- |                              |                      |
|------------------------------|----------------------|
| a. Mains Adapter             | j. SINGLE button     |
| b. Keys                      | k. Hand-held Wand    |
| c. DC Power Socket           | l. Laser Button      |
| d. Mute Button               | m. Orange Light      |
| e. Keyswitch                 | n. White Light       |
| f. Power Level Indicator     | o. Touch Pins        |
| g. Power Increase Button (+) | p. Laser Aperture    |
| h. Power Decrease Button (-) | q. Indentation Guide |
| i. SCAN button               | r. Front Face        |

## treatment area


### "SCAN" mode

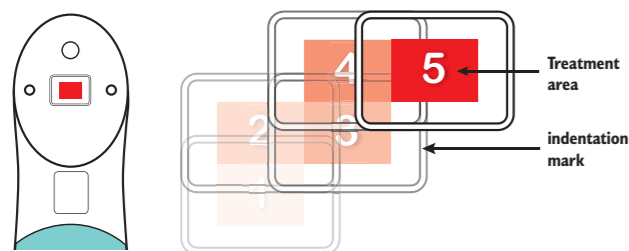

As a guide to the treatment area, a small indentation mark will be left by pressing the Wand down on the skin. In SCAN mode, the laser covers an area of 7mm x 5mm.

Overlapping the Indentation mark by a third when moving the Wand in each direction will ensure effective coverage.

### "SINGLE" mode

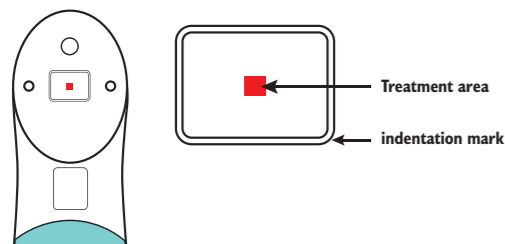

SINGLE mode is for treating individual hairs. The treatment area is in the centre of the indentation guide and measures 2mm x 2mm.

## laser hair treatment

### patch test

**IMPORTANT:** 24 hours before treatment it is very important to do a patch test.

The patch test will confirm the maximum power level for your skin type. Each area to be treated must be checked as skin colour can vary. Please refer to the 'PATCH TEST' section in the Information Guide and on the DVD for instructions.


A slight redness after treatment may occur. Any localised redness after treatment is caused by the destruction of the hair follicles and is quite normal. However, if any redness remains after 24 hours this would indicate the power level selected was too high for your skin type.

### getting started


Please ensure you read the Information Guide and understand the 'Warnings' and 'Product Safety' sections before proceeding.

For best results, it is recommended you use the SCAN mode. This does not require the individual targeting of a hair because it treats an entire area.

- Best results are achieved on hairs that have been trimmed to 1mm–2mm in length.
- Remove all jewellery in the vicinity of the treatment area.
- Please ensure you have patch tested the area to be treated 24 hours prior to treatment.

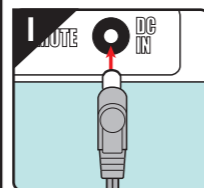

1. Connect Mains Adapter into the DC Power Socket on the back of the unit and plug into the Mains. The Keyswitch will light up orange.

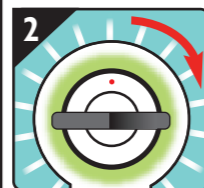


2. Unlock with Key to turn on. The light around the Keyswitch will flash green.




3. Enter the "SECURITY BUTTON SEQUENCE" within 60 seconds of turning the Keyswitch ON.

You will find the Security Button Sequence in the Information Guide.

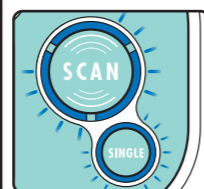

4. The green light around the Keyswitch will fully illuminate. The lights around the SCAN and SINGLE buttons will flash blue ready for selection. The bottom Power Level Indicator light will illuminate.

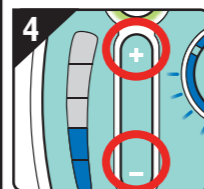


5. Use the +/- buttons to select the correct power level for your skin type. Refer to the 'Patch Test' section of the Information Guide for skin types.

Each time the power setting is adjusted, you will need to reselect SCAN or SINGLE mode to continue.




6. If no buttons are pressed after 20 minutes the Unit will power off completely. To restart turn key off and back on again and start at Step 3.

### 'SCAN' hair treatment

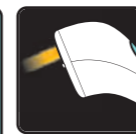

5. Press the SCAN button once. The blue light around the SCAN button will fully illuminate. The orange light on the Wand will illuminate.

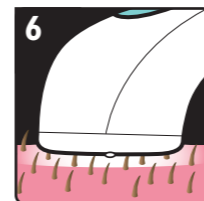


6. Place the Wand on the skin over the area you wish to treat, ensuring the Touch Pins are in contact with the skin.

Refer to 'treatment area' section, SCAN mode.

The orange light will turn OFF and the white 'ready' light will turn ON.




The blue lights around the SCAN button will rotate to indicate treatment is ready to start.

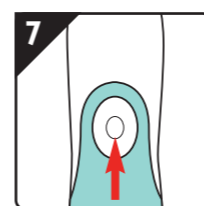


7. Press and HOLD the Laser Button on the top of the wand to start laser treatment.

The white light will flash on and off to indicate laser treatment is in progress and the unit will beep repeatedly (this can be muted). You will not be able to see the laser.

When the laser treatment ceases, the orange light will come back on. You can now move onto the next area of hairs to be treated.

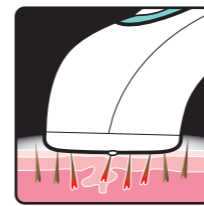

You may feel a 'pin prick' sensation during treatment. Any localised redness after treatment is caused by the destruction of the hair follicles and is quite normal. However, if any redness remains after 24 hours this would indicate the power level selected was too high for your skin type. Please refer to the 'Patch Test' section in the Information Guide. Also note, any longer lasting redness will also disappear.

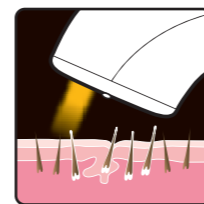

When moving the Wand to the next area of hairs to be treated the Touch Pins must be removed briefly from the skin.

Refer to 'treatment area' section, SCAN mode.



Do NOT move the Wand during treatment. You may notice the Wand becomes warm during use this is quite normal.



In SCAN mode, the laser will run for 10 seconds on the lowest power level, and up to 50 seconds on the highest level. In SINGLE mode the laser will run for 2 seconds on the lowest power level, and up to 10 seconds on the highest power level.



If contact is lost between the Touch Pins and skin during treatment, the orange light will come on. Treatment will resume when contact is made again.



Sometimes you can smell faint burning or even see a very, small line of smoke as the hair is treated.



Successfully treated hairs will fall out within 7–14 days (light rubbing will encourage this). See 'Hair Growth Cycle' section in the Information Guide.




8. Turn the Keyswitch anticlockwise to lock the unit. Remove the Key to a safe place away from children. The Keyswitch will illuminate orange until the adapter is unplugged.



If you experience unexpected behaviour of this product then it is possible some external electromagnetic interference is responsible e.g. electric motors/loudspeakers etc. You should remove this source and restart the product to restore full operation. At no time will the laser fire unintentionally.

### 'SINGLE' hair treatment

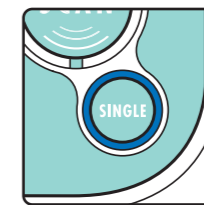


You can opt to treat single hairs using the SINGLE mode.

Press the 'SINGLE' button once. The SINGLE button will illuminate blue.

The orange light on the Wand will illuminate.

Follow the same procedure for SCAN treatment from Step 6, but this time pressing the SINGLE button.

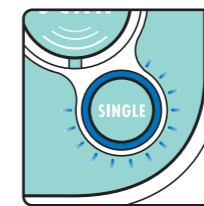


In SINGLE mode, the blue light around the SINGLE button will flash.

Treatment time is between 2–10 seconds in SINGLE mode depending on the power level selected.

Refer to 'aiming and positioning' section, SINGLE mode.

### additional information

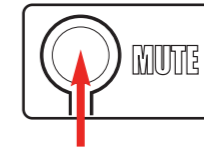

Press the MODE button on the back of the Unit to turn off the Beep noises. Press again to turn back on.


This doesn't turn off warning beeps.

### resetting after 'lockout'



If you experience a 'lockout' as indicated by repeated flashing of any of the blue power lights, simply reset by turning the unit OFF using the key, waiting 5 seconds and turning the unit on again.

### cleaning laser aperture

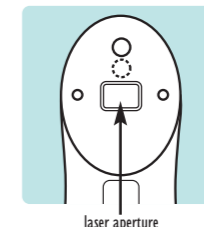

Periodically check whether there are any smoke deposits on the Laser Aperture. A build up of deposits will affect the efficiency of the laser. If you notice any deposits, switch the unit OFF and unplug the mains adapter. Then simply clean with a lint free cloth or cotton bud, ensuring no fibres are left in the Laser Aperture.