

SALON LASER

HAIR REMOVAL SYSTEM



808nm wavelength

Wearing laser safety goggles is unnecessary as the product has been designed with inbuilt safety features.

LASER SAFE TECHNOLOGY

PATENTS PENDING



CLASS 1 PRODUCT CONTAINING CLASS 3B LASER

1-IN-LAHR-UK/3.0



This product **MUST** be used responsibly.

Please read **ALL** instructions including the User Manual **CAREFULLY** before use.



Refer to your User Manual for:

- SECURITY BUTTON CODE
- INFORMATION ON HAIR, LASERS, SKIN TYPES AND TREATMENTS
- AFTERCARE ADVICE
- HANDY TIPS
- UNIT MAINTENANCE
- TROUBLESHOOTING
- QUESTIONS & ANSWERS
- CUSTOMER SERVICE

INSTRUCTIONS

WARNINGS

Not suitable for Afro-Caribbean or black skin complexions. See "Suitable Skin Types" in User Manual.

Laser hair removal should not cause discomfort. If you feel any discomfort or soreness the setting may be too high or your skin tone (colour) is not suitable for use with this product.

Inbuilt safety features protect the eye being exposed to the laser, however, never aim the Wand at or near eyes.

Do NOT use near flammable materials.

Keep the product out of reach of children and pets.

Reflection Hazard: Please remove all jewellery in the vicinity of the treatment area before use.

When not in use, store the key responsibly in a safe place away from children.

Do not treat the same hair repeatedly in any one session.

Do not use after drinking alcohol or if you are feeling unwell or fatigued.

Do not treat the following areas: tattoos, permanent make-up, pimples, warts, hair on black spots or moles, dark freckles, skin disease, birth marks, burns, lips, nipples, genital or anal areas, eczema, eyelashes, eyebrows, near the eyes, nasal or ear hair.

CAUTIONS


Do NOT dismantle or repair. If you have any problems using the product, contact our Customer Care Department.

Misuse of the unit can cause harm.

Do NOT override any safety features or open the unit

Only use this product with the Mains Adapter supplied.

COMPONENTS

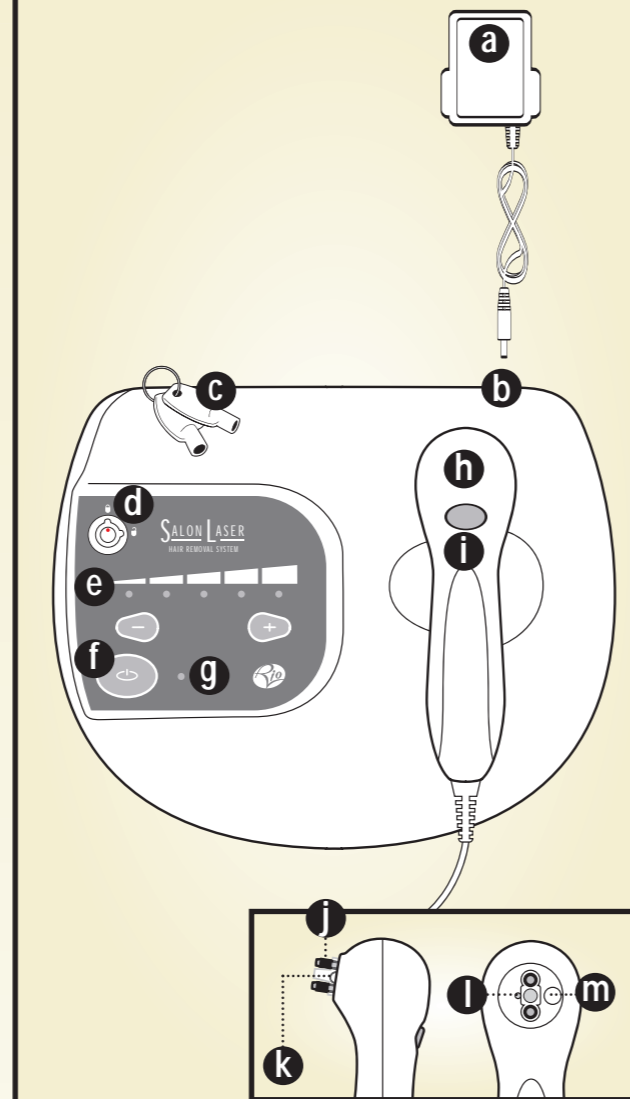

- | | |
|---------------------------------|--------------------------------|
| a. Mains Adapter | h. Hand-held Wand |
| b. Power Socket | i. Laser Button |
| c. Keys | j. Laser Shroud and Touch Pins |
| d. Keypad | k. Laser Aperture |
| e. Power Level (5 green lights) | l. White 'ready' Light |
| f. Enable Button | m. Amber Light |
| g. Enable Light (red) | |

PATCH TEST


IMPORTANT: 24 hours before treatment it is very important to do a patch test.

The patch test will confirm the maximum power level for your skin type. Each area to be treated must be checked as skin colour can vary.


Please refer to the 'PATCH TEST' section in the user manual for instructions.

TREATMENT


- Best results are achieved on hairs that have been shaved a day or two before treatment or trimmed to 1mm-2mm in length.
- Remove all jewellery in the vicinity of the treatment area.
- Please ensure you have patch tested the area to be treated 24 hours prior to treatment.

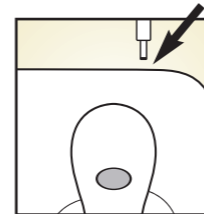


1. Plug Mains Adapter into socket on the Unit and into the Mains.

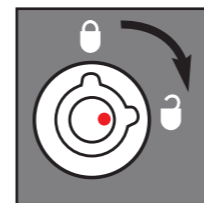


2. Use the Key to unlock the Keypad. The left Power Level light will flash green.

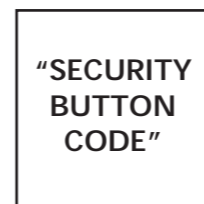


3. Enter "SECURITY BUTTON CODE" within 60 seconds. You will find the Security Button Code in the front of the User Manual.

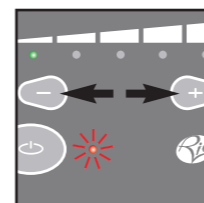

4. The + and - buttons increase or reduce power levels. Refer to 'Laser Power Levels' in the User Manual for the maximum power setting for your skin type.

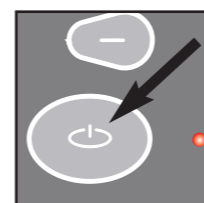


- If you feel any discomfort reduce the power level.
- If the power level is adjusted the red Enable light will flash.
5. Press the Enable Button. The red Enable light will stay on and the amber light on the Wand will illuminate.

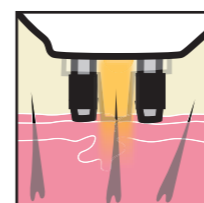


6. **Dampen** the area to be treated with water. Place the Touch Pins either side of the hair you wish to treat.

Important: The unit will not work if the skin is too dry.


TREATMENT

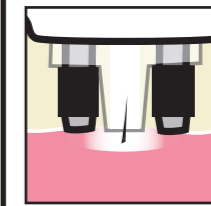


7. The amber light will go OFF and the white 'ready' light will illuminate between the Touch Pins. **Important: If the amber light does not go out, apply more water or use a water based moisturiser.**

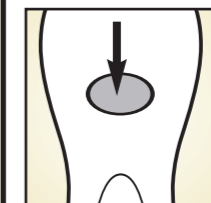


8. Press and HOLD the Laser Button to start the treatment. **This must be aimed at the hair ROOT.** You will hear a series of high pitch beeps to indicate the laser is operating. The laser appears as a small red light.

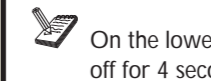


- On the lower power levels, the laser will pulse on and off for 4 seconds.

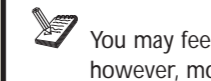

- You may feel a 'pin prick' sensation during treatment, however, most people do not feel anything.

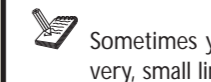

- Sometimes you can smell faint burning or even see a very, small line of smoke as the hair is treated.

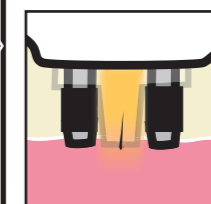


9. When the red laser light treatment ceases the amber light will come back on. You can now move onto the next hair.



- The Touch Pins must be removed briefly from the skin between treating each hair.



- Successfully treated hairs will fall out within a **few days** (light rubbing will encourage this).

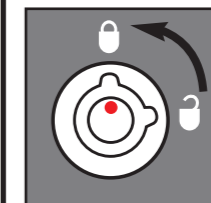

11. After use, turn the Keypad to lock the unit. Remove the Key to a safe place away from children.



- If no buttons are pressed within 60 seconds the unit will enter a 'resting mode'. After 20 minutes the Unit will power off completely.

If you experience unexpected behaviour of this product then it is possible some external electromagnetic interference is responsible. You should remove this source and restart the product to restore full operation. At no time will the laser fire unintentionally.